CheckToP test kit



CheckToP is a urine test kit you can use at home to check if a pregnancy has ended after an early medical abortion.

- The test checks your urine (pee) for a hormone called hCG.
- Your body makes hCG when you're pregnant.
- The amount of hCG in your body decreases when you're not pregnant anymore.

Note: hCG = human chorionic gonadotropin



When to do the CheckToP test


In an early medical abortion, you're given two medicines to end a pregnancy – mifepristone and misoprostol.



3 to 4 weeks after taking these medicines, the CheckToP test is used to confirm the pregnancy has ended.

Timing is important


- **Don't test too early (before** 21 days). Your hCG levels could still be high.
- **Don't test too late** (after 28 days). You'll delay extra care and support if the medicine hasn't worked.

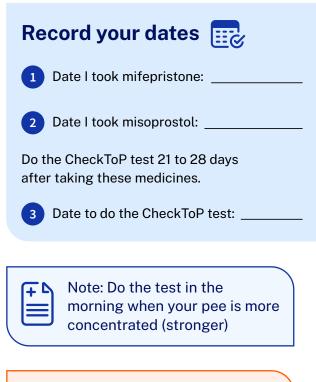

When done correctly, the CheckToP test is very reliable (99%).

How to do the test

1. Take the cap off the stick.

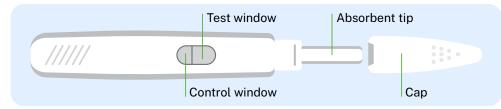

- 2. Pee on the absorbent tip for 5 to 10 seconds, **OR** pee in a clean, dry container and dip the tip in for 5 to 10 seconds.
- **3.** It's important you don't get any pee on the windows this could affect the result.
- **4.** Put the cap back on and lay the stick on a dry flat surface with the windows facing up.
- 5. Wait between 5 to 10 minutes to read the test results.




Important: Do not read the result before 5 minutes or after 10 minutes – It may affect the results.

For more information about the CheckToP test, visit healthify.nz/checktop

Reading your results

Note: The CheckToP test stick has 2 windows:

- a test window (larger window)
- a control window (smaller window).

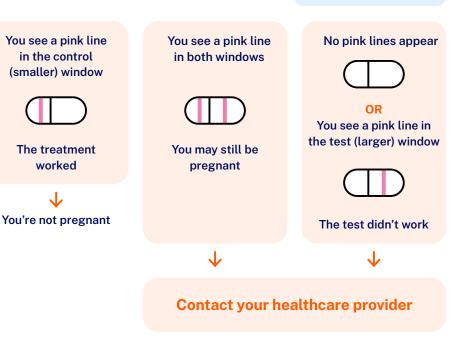

Signs and symptoms you might still be pregnant:

- little or no bleeding after taking misoprostol
- nausea (feeling sick)
- sore breasts
- tummy is growing
- it's been 4 to 6 weeks and you haven't had a period.



Test window

Control window